

INTRODUCTION

Micrometeorites are tiny extraterrestrial particles, which survived atmospheric entry. They fall to Earth at a rate of 40,000 tons annually, and can be retrieved e.g. from the Antarctic. They represent a valuable source of information on the chemical evolution of the Solar System. The microscopic size (50-2000 μm) of micrometeorites requires the use of novel *in situ* mapping techniques of elemental analysis: with high lateral resolution, preferably non-destructive, with multi-element capabilities, and capability to provide quantitative data. Hyphenated with a low-dispersion laser ablation (LA) system, time of flight ICP-MS (ICP-ToF-MS) allows for 2D mapping with a laser repetition rate of several hundred Hz, with each laser shot recorded as a single pixel, and rapid quasi-simultaneous acquisition of almost the entire periodic table.

AIMS AND GOALS

- Evaluation of the applicability of **LA-ICP-ToF-MS method for 2D element mapping** analysis of 50-2000 μm $\varnothing$ cosmic spherules (CS)
- Development of a **calibration procedure** for fully quantitative mapping (using GeoReM consensus values)
- **Validation using glass reference materials** of natural composition

SAMPLES

- **Melted micrometeorites (cosmic spherules, CS)** larger than $\sim 200\mu\text{m}$ (glassy, barred olivine, cryptocrystalline, and Ca-Al-Ti rich types)
- Collected from sedimentary traps near the Widerøefjellet (altitude $\sim 2750\text{m}$), Sør Rondane Mountains, **East Antarctica**
- Casted into epoxy mounts and polished with 4000 grit diamond paste

LA-ICP-ToF-MS METHOD VALIDATION

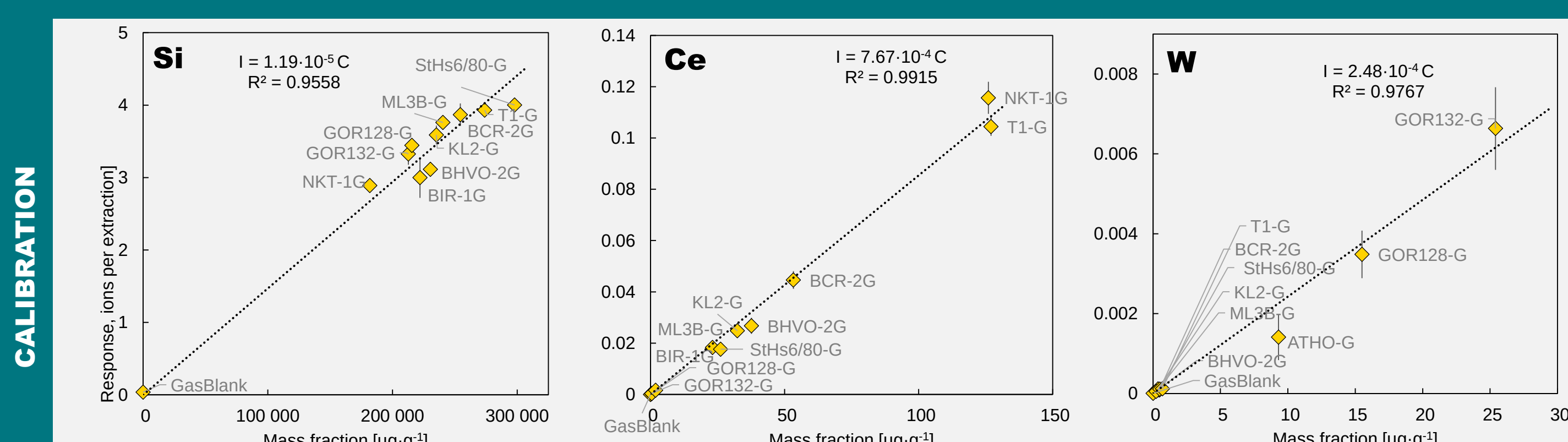


Fig. 1. Calibration lines for the selected elements using LA-ICP-ToF-MS mapping

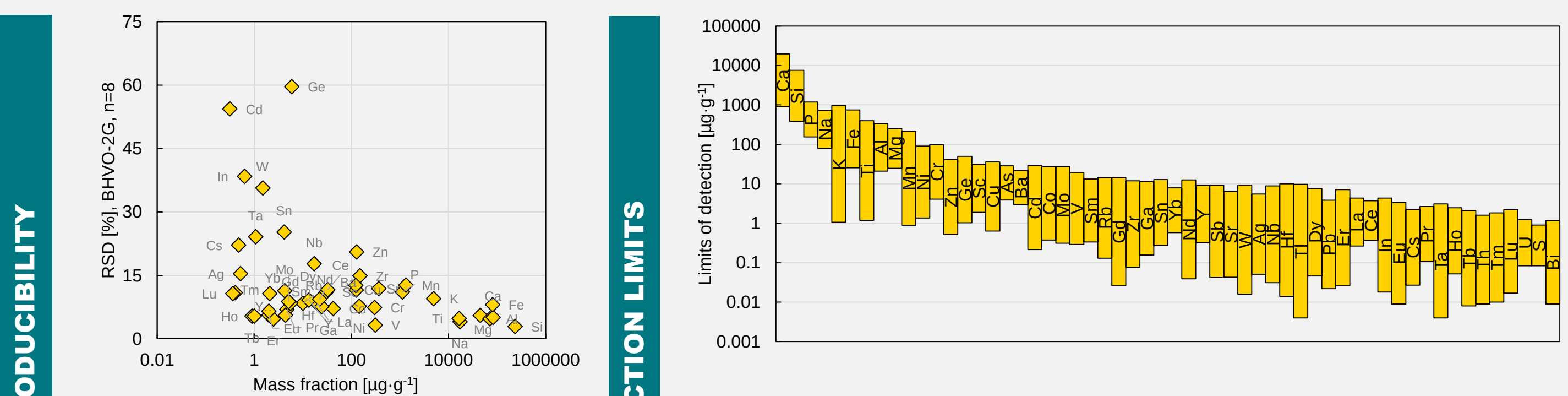


Fig. 2. Intermediate precision of the LA-ICP-ToF-MS mapping (results of USGS BHVO-2G over a period of ~ 1 year)

Fig. 3. Limits of detection of the LA-ICP-ToF-MS mapping for integrated areas according to the equation of Pettke *et al.*, 2012 (3x3 μm min-max range over a period of ~ 1 year)

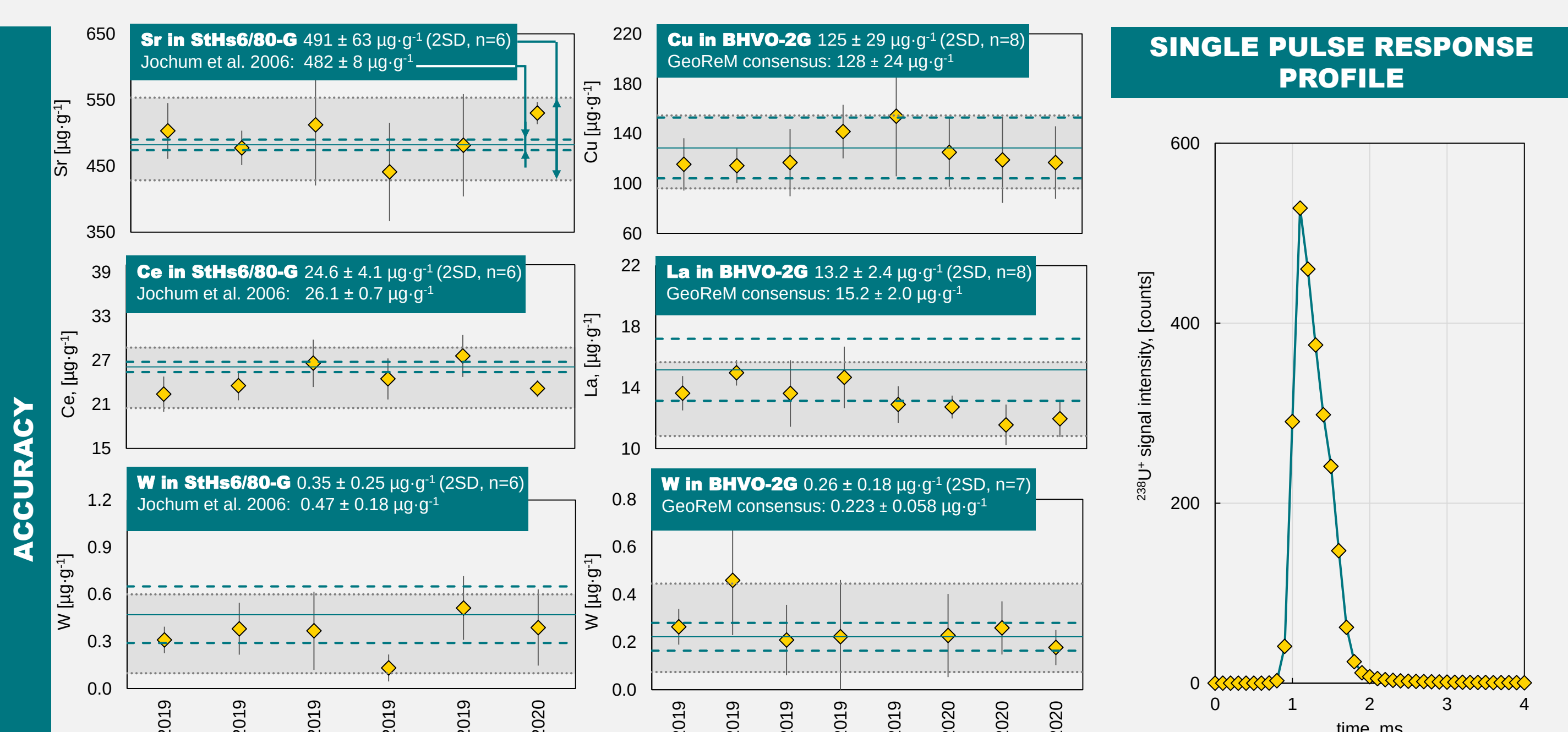


Fig. 4. LA-ICP-ToF-MS mapping of glass RMs shows no bias and decent intermediate precision

SINGLE PULSE RESPONSE PROFILE

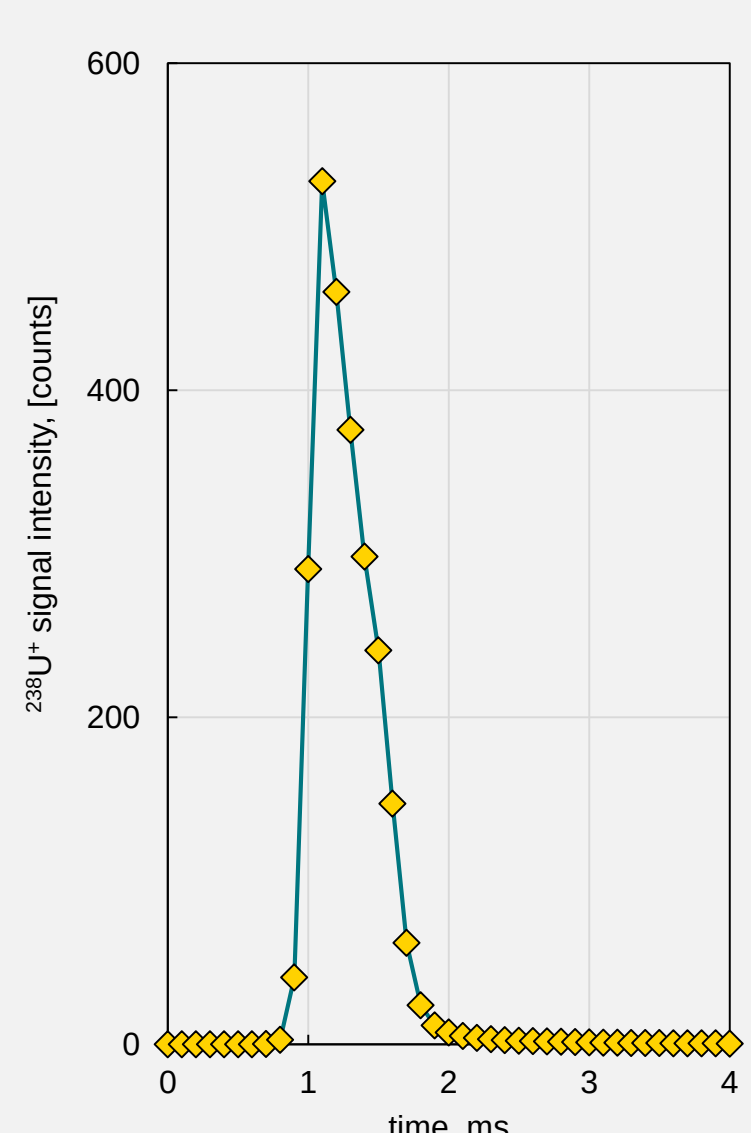


Fig. 5. Single pulse response profile of the LA-ICP-ToF-MS setup


ACKNOWLEDGEMENTS

S.M.Ch. acknowledges his postdoctoral fellowship from the Flemish Research Foundation (FWO, "Excellence of Science" project ET-HoME - ID 30442502) and the support from the Montanuniversität Leoben for attending the EWCPs. T.V.A. acknowledges the FWO for his postdoctoral research fellowship (FWO.3E0.2022.0048.01). Collection of micrometeorites was made possible by the Belgian Science Policy through financial support for various projects (BELAM, Amundsen and BAMM). F.V. acknowledges Teledyne Photon Machines for logistic and financial support.

METHODS

- **icpTOF2R (TOFWERK) Time-of-Flight ICP-MS** unit. Quasi-simultaneous detection (ToF) of a m/z range 23 - 238 ($^{23}\text{Na}^+$ to $^{238}\text{U}^+$), base integration time 48 μs , $m/z = 40$ ($^{40}\text{Ar}^+$) notch-filtered
- **IRIDIA (Teledyne Photon Machines)** 193 nm ArF* laser ablation system, Cobalt™ ablation chamber, ARIS aerosol rapid introduction system. Mapping speed: 0.3 - 30 $\text{mm}^2 \cdot \text{h}^{-1}$, spot size 1x1 to 5x5 μm , laser repetition rate 100-300 Hz. Single pulse response duration ~ 1 ms, fluence: ~ 2 $\text{J} \cdot \text{cm}^{-2}$
- **HDIP** processing software for visualization, image calibration and segmentation
- Multi-point calibration using USGS and MPI-DING glass reference materials (GSE-1G, GSD-1G, BCR-2G, NKT-1G, BIR-1G, BHVO-2G, GOR132-G, GOR128-G, ML3B-G, ATHO-G, StHs6/80-G, KL2-G, T1-G)
- **GeoReM consensus values** are relied onto where certified data is missing
- Normalization of the oxides to 100 % to account for ablation yield variations

MAPS OF THE ELEMENT MASS FRACTIONS

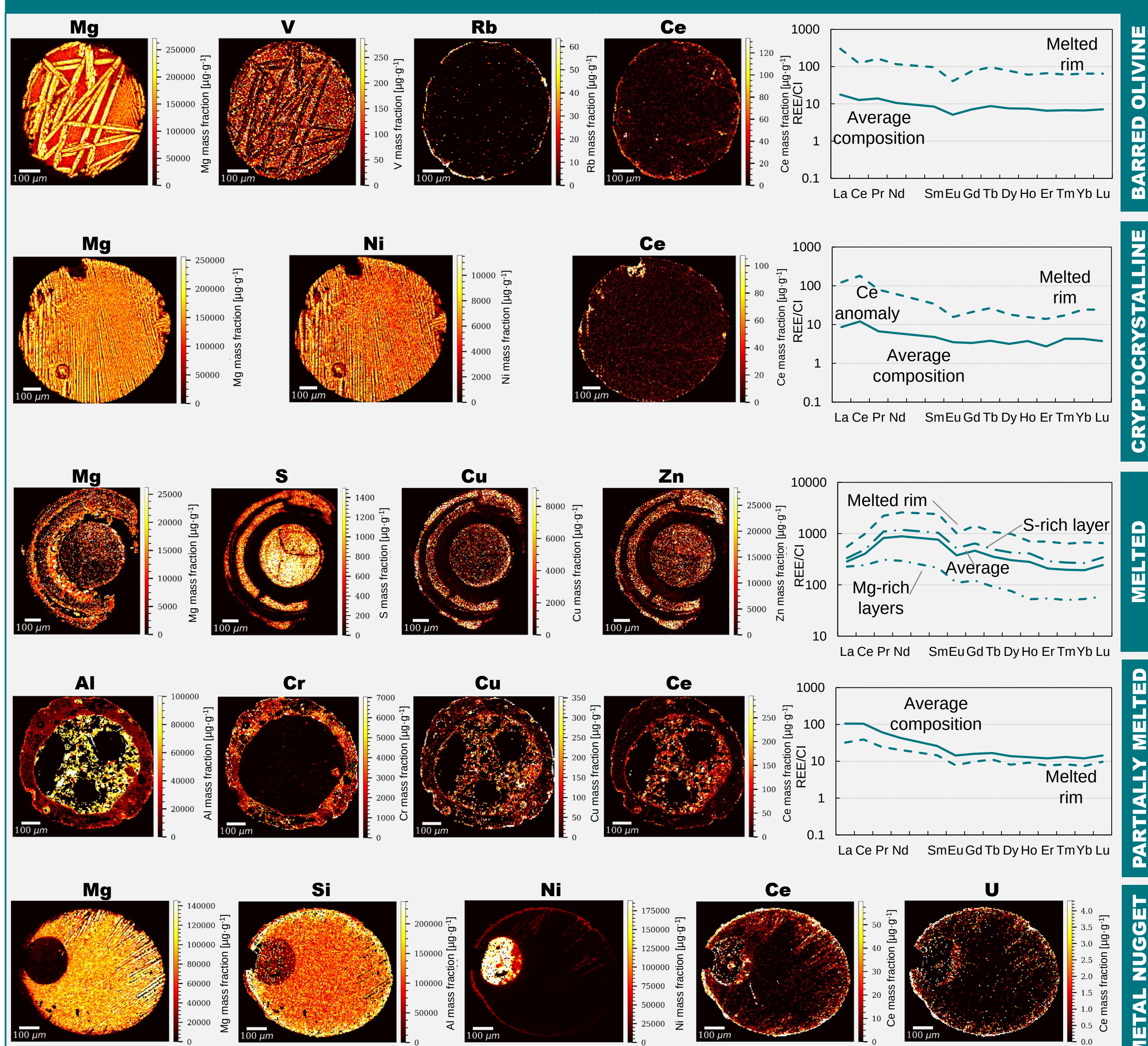

CONCLUSIONS AND OUTLOOK

- Element range from **Na to U** is covered **for each laser shot** (at 100-300 Hz frequency each shot corresponds to a single pixel 1x1 to 5x5 μm)
- **Fully quantitative element maps**
 - Multi-point calibration
 - Matrix matched glass RMs
 - Reliability of the 100% oxide normalization is better than that of an internal standard
- **Calibrated element data has good precision & no systematic bias**
 - Absence of bias is verified using glass RMs of natural rock composition
 - **Intermediate precision ~ 5 -15% (SD)** for most elements $> 10 \mu\text{g} \cdot \text{g}^{-1}$, and ~ 5 -30% (SD) for the elements $< 10 \mu\text{g} \cdot \text{g}^{-1}$
 - Data for some elements (e.g. Cd, Ge, W, In, Sn, As) is only semi-quantitative due to low intensities, poor data for RMs, non-homogenous RMs (?)
 - **Detection limits of ~ 0.1 -10 $\mu\text{g} \cdot \text{g}^{-1}$** for integrated areas for element maps with 3 μm spot size (~ 100 -10000 $\mu\text{g} \cdot \text{g}^{-1}$ for major elements)
- **The LA-ICP-ToF-MS maps compliment petrography observations**
 - Melts of different composition can be recognized, segmented and quantified
 - The maps and element patterns confirm the extraterrestrial nature of the CS and inform on the processes of atmospheric melting and terrestrial residence.

